

BladeSINGLE

One | Blade Remover



A *simply safer* way of removing blades in operating rooms.

Remove, contain and dispose of used blades without the unacceptable sharps injury risk of two-handed removal techniques involving fingers, forceps or re-sheathing.

The surgeons have their preferred steel handle and the nurses are protected.



The essentials of team safety

A simply safer way to stay focused on your job.

Clicksmart Blades**SINGLE** provides a *simply safer way* to remove and contain used scalpel blades in sterile operating environments. This world-first single-handed scalpel blade remover prevents the dangerous injuries commonly caused by removing blades with fingers, forceps or re-sheathing techniques.

Continue to use the best tool for your job.



Most surgeons prefer the balance, weight and ergonomics of a reusable steel scalpel but removing used blades manually compromises staff safety.

Flimsy disposable plastic scalpels are not the answer because they are harder to use, too light to balance, too cumbersome to use, obstruct the line of sight and get entangled with the tissues. It's like eating your dinner with a plastic knife and fork.

Make compliance a practical reality.

Research shows that when a single-handed scalpel blade remover is combined with a hands free passing technique (HFPT) it is up to five times safer than a safety scalpel.⁴

Blades**SINGLE** meets these essential safety standards:

- OSHA Bloodborne Pathogen Standard "1910.1030"
- Australian/New Zealand Standard "AS/NZS 3825:1998"
- CE compliant

This essential safety solution makes compliance a practical reality and dramatically reduces injury costs. Healthcare workers in the United States alone report between 600,000 and 1 million sharps injuries per year. Many more go unreported.⁵ Occupational health and safety regulations worldwide are responding.


Blades**SINGLE** ensures optimal safety based on speed and convenience for the end-user, ergonomic simplicity and passive, automatic and single-handed safety principles.

Minimise risks of injury and infection


Sharps injuries are a serious occupational hazard in healthcare settings. Scalpel cuts are not only more severe than needle stick injuries, they are also more likely to happen. The single-handed technique of removing blades minimizes your risk.

- 8-11% of sharps injuries are caused by scalpel blades.^{1,2}
- There are 662 scalpel injuries per 100,000 scalpel blades used, compared to 3.2 syringe and loose needle incidences per 100,000 disposable syringes and loose needles used.³ Scalpel injury rates are 200 times higher than sharps injuries involving syringes.


These alarmingly common scalpel injuries cause pain, emotional distress, potentially fatal infections, and damage to digital nerves, arteries or tendons requiring microsurgery and extensive rehabilitation. Blades**SINGLE** dramatically reduces the risk of scalpel blade injuries in operating theatres.

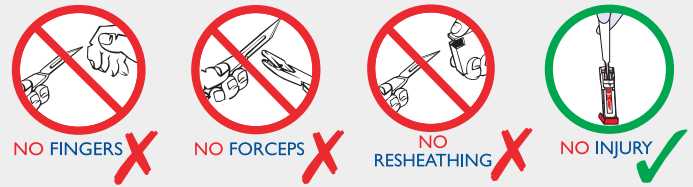

Devices Involved in Percutaneous Injuries¹

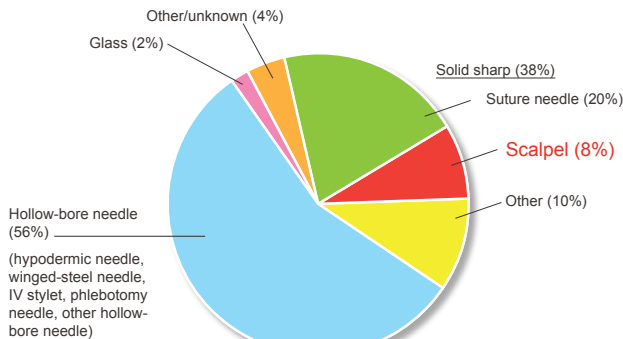

Incidences per 100,000 devices used³

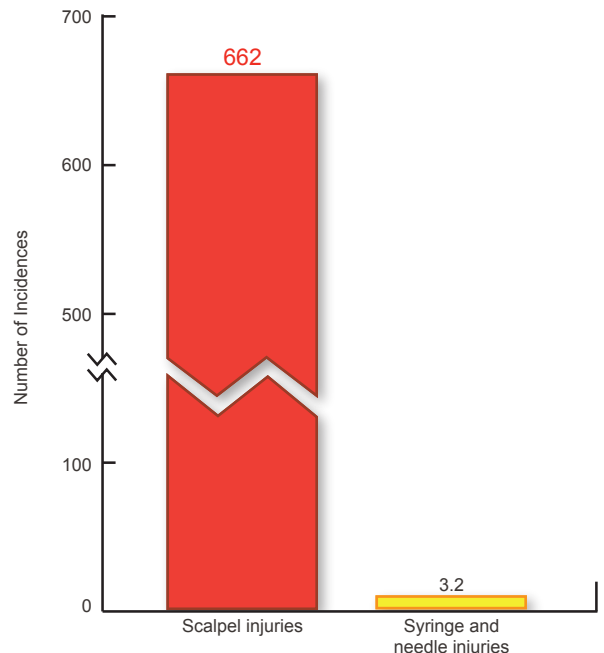

1. "CDC (Centre for Disease Control and Prevention) Sharps Injury Prevention Workbook 2008." from <http://www.cdc.gov/SharpsSafety/index.html>
 2. Jagger, J., M. Bentley, et al. (1998) "A study of patterns and prevention of blood exposures in OR personnel." AORN Journal 67(5): 979-81, 983-4, 986-7 passim.
 3. Eisenstein, H. C. and D. A. Smith (1992). "Epidemiology of reported sharps injuries in a tertiary care hospital." J Hosp Infect 20(4): 271-80.
 4. Fuentes et al. (2008) "Scalpel Safety: Modeling the effectiveness of different safety devices' ability to reduce scalpel blade injuries." The International Journal of Risk & Safety in Medicine 20(1-2):83-89.
 5. Matson, K. (2000) "States begin passing sharps and needle-stick legislation to protect health care workers." AORN Journal 72(4): 699-703, 705-7.

BladesINGLE


One | Blade Remover

Remove and contain blades with one hand - push and click

With the single-use BladesINGLE you can simply 'push and click' to change and contain blades during operating without risking the other hand. Using double-sided tape, the red BladesINGLE base grips on the drapes and avoids slipping, making single-hand removal of blades possible. The other hand is out of harm's way until the blade is safely contained in the remover.

The BladesINGLE simple 'push and click' blade removal solution for steel scalpels:

- **IMPROVES PATIENT SAFETY** by enabling surgeons to use ergonomic precision steel scalpels.
- **IMPROVES STAFF SAFETY** by ensuring the used scalpel blade is removed with a safe single-handed technique and then safely contained until final disposal.

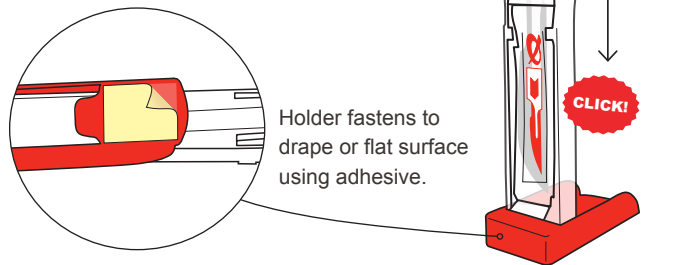

Contain blades immediately after use

Contaminated scalpels lying around on the mayo table increases the risk of exposure to sharps injury and potential bloodborne pathogens. The sterile BladesINGLE removes and contains contaminated blades for safe disposal immediately after use, eliminating lurking sharps hazards.

Ensure team safety and satisfaction

According to EPINet data in 2001, percutaneous injury rates show that nurses are most at risk (43.6%) of sustaining sharps injuries.⁹

With BladesINGLE, surgeons use their preferred traditional steel scalpels while nurses are protected.

Sterile device

BladesINGLE can be included in sterile packs (e.g. custom kits or theatre packs), ensuring that the blade remover is immediately available for use when required. BladesINGLE can be sterilised by ETO or gamma radiation.

Scalpel Safety vs. Safety Scalpel

EPINet (Exposure Prevention Information Network) data published in 2003 by Perry et al. noted that in the year 2000-2001, there were 91 injuries caused by reusable scalpels and 42 injuries caused by disposable scalpels (described elsewhere in the text as safety scalpels).⁶

In the year 2000-2001, there was only a 10% market share of safety scalpels while 90% of the population were still using the traditional handle.

The conclusion would therefore be that the relative incidence of injuries is four times higher with the use of safety scalpels than it is for traditional reusable scalpel handles.

CDC, in a landmark study on active and passive devices, found that up to 90% of active devices had not been activated.⁷ Safety scalpels are active devices, and this could explain the higher relative incidence of injury caused by safety scalpels.

Fuentes et al. reviewed the largest data set of scalpel injuries published at that time and found that the combination of a single handed scalpel blade remover when used with a hands free passing technique was as safe and up to five times safer than a safety scalpel.⁴ This is **scalpel safety**.

OSHA updates


OSHA in an update on December 22nd, 2005 stated:⁸

"... if the use of a particular engineering control, in this case a safety scalpel, compromises patient safety, its use would not be considered feasible... blade removal must be accomplished through the use of a mechanical device or one handed technique... The use of a single-handed scalpel blade remover meets these criteria..."

Counting of blades

BladesINGLE has been specifically designed to enable easy and accurate counting of blades after use in theatre.

Fits seamlessly into your standard operating room procedures

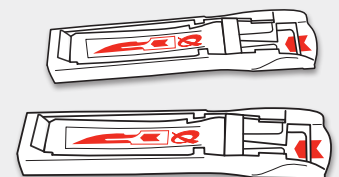

6. Perry, J. and J. Jagger (2003). "Scalpel Blades: Reducing Injury Risk." *Advances in Exposure Prevention* 6(4): 37-48.

7. Alvarado-Ramy, F., E. M. Beltrami, et al. (2003). "A comprehensive approach to percutaneous injury prevention during phlebotomy: results of a multicenter study, 1993-1995." *Infect Control Hosp Epidemiol* 24(2): 97-104.

8. OSHA Standard Interpretations. Standard Number: 1910.1030; 1910.1030(c)(1)(iv); 1910.1030(d)(2). *12/22/2005 - Use of passing trays and single-handed scalpel blade remover in a surgical setting.*

9. Uniform Needlestick and Sharp-Object Injury Report U.S. EPINet Network, 2001, 58 Healthcare facilities.

Specifications

Product code	QSBS01	
GMDN code	46236	
Blade capacity	1 blade - Each Blades ^{SINGLE} cartridge removes 1 (one) blade.	
Compatible Blades & Handles*	Handles: #3, 3L, 4, 4L, 5, 7, 9, baron, cylindrical. Blades: #10, 10a, 11, 11P, 12, 14, 15, 15a, 15c, 16, 19, 20, 21, 22, 22a, 23, 24, 25, 25a, 27, 36, 6, 9, E11, Sabre D/15.	
Incompatible Blades & Handles	Beaver handle, Disposable handles	
Size (per unit)	Folded: Length 9.3cm Width 2.2cm Thickness 0.9cm	Standing: Length 7.3cm Width 2.2cm Height 9.3cm
Weight (per unit)	< 10g	

* Most brands will work with Clicksmart Blades^{SINGLE}

Material	Borealis polypropylene (PP)
Latex	No
DEHP	No
Sterile	Yes
Sterilisation method	Gamma radiation or ETO
Re-usable	No, each unit removes 1 (one) blade before disposal.
Disposal method	After the blade count has been performed, the whole Blades ^{SINGLE} is disposed of as medical waste in a sharps container.
Regulatory approvals	Australia - Listed under TGA ARTG NO: 134623 USA - FDA listed Europe - Class I Sterile medical device under Council Directive 93/42/EEC and CE marked


Instructions for Use

Blades^{SINGLE}


One | Blade Remover

PREPARE USING TWO HANDS:


1. Peel to reveal adhesive.

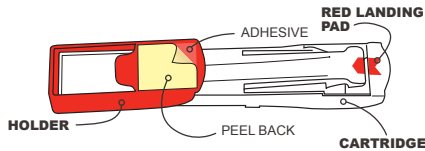

2. Take cartridge out of holder.

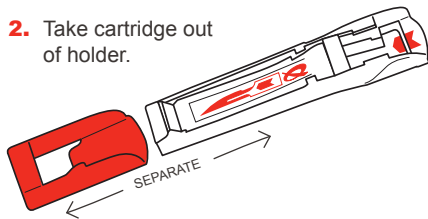

3. Fasten holder to drape or flat surface.

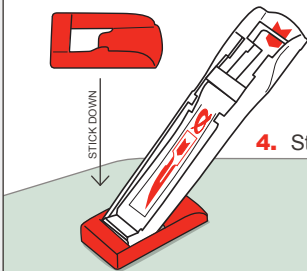


4. Stand in holder.


READY FOR USE

ACTIVATE SINGLE HANDED:


5. Using one hand, slide tip of scalpel into cartridge.

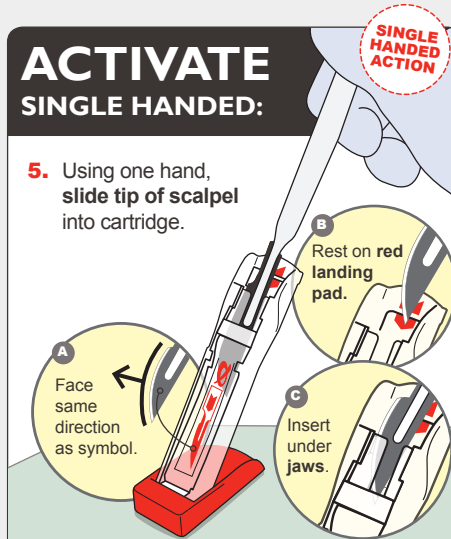

6. **TILT UP** then push down.

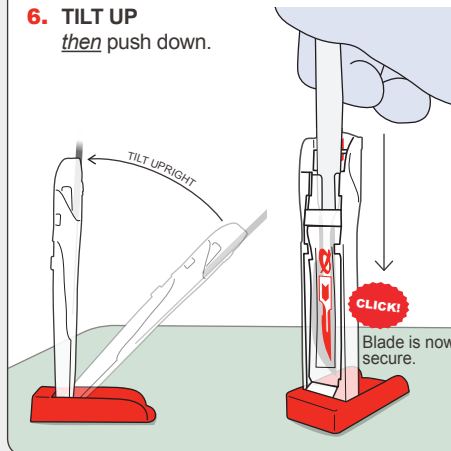


7. Remove cartridge from holder.

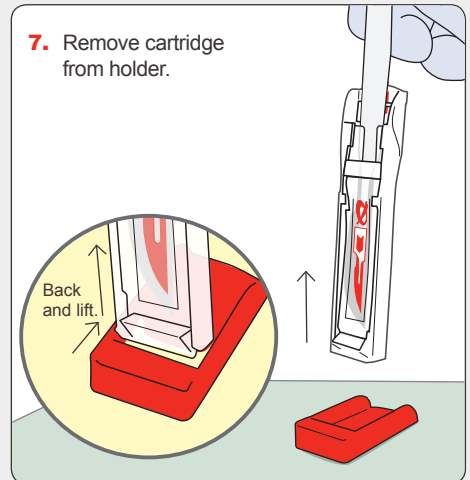

SAFE TO USE TWO HANDS:

8. Pull out scalpel handle.

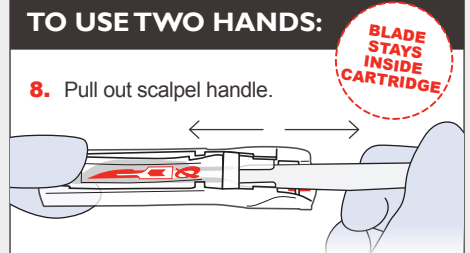

9. Ready for counting and later disposal in sharps bin.



QLICKSMART Pty Ltd.

+61-7-3844-1182 T

+61-7-3844-1183 F

mjs@clicksmart.com

www.clicksmart.com

 **QLICKSMART**[®]
Simply safer